


BioProfile® FLEX
The Power of One

One Automated Modular Analyzer for Fast Comprehensive Cell Culture Analysis

WITH
NEW OPC
CONNECTIVITY
SUITE

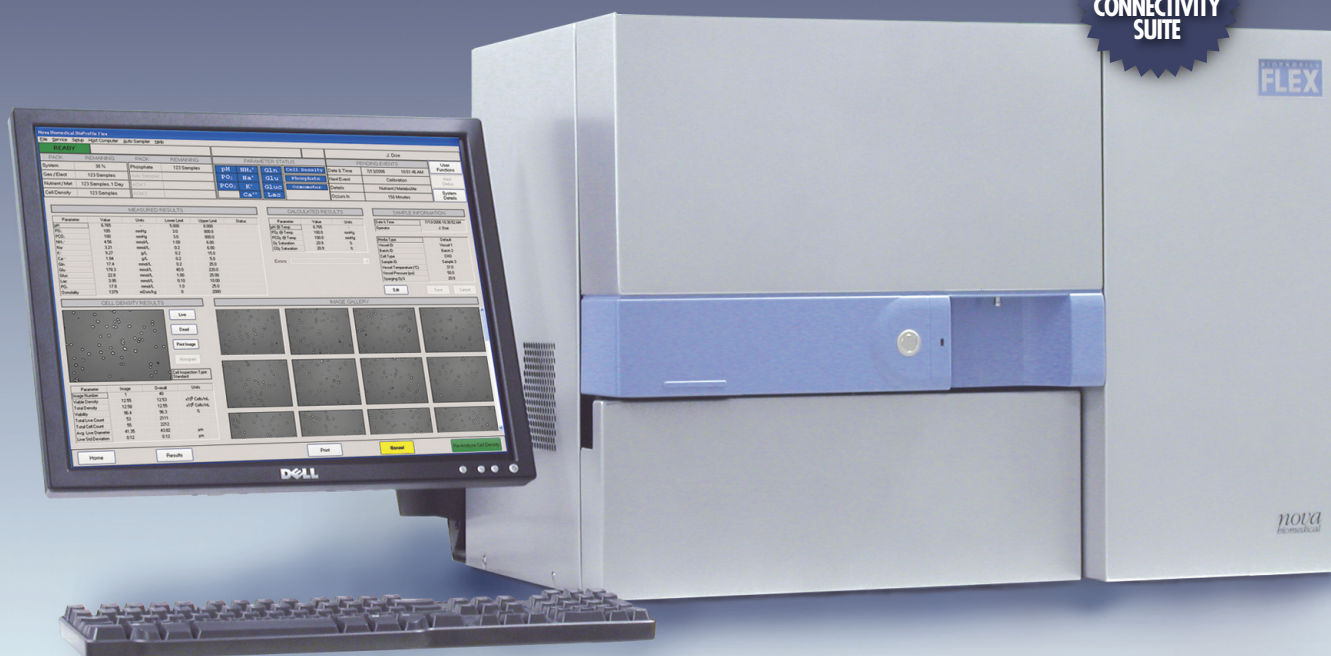

The Power of One

BioProfile FLEX can reduce cell culture analysis time, labor, and operating costs by consolidating multiple analyzers into a single, easy to use workstation.

- One compact workstation versus 4 or 5 separate analyzers
- One small, 1 mL sample conserves cell culture mass and end product
- One fast, 6 minute analysis saves at least one hour over multiple instrument analyses
- One integrated data source simplifies data collection, analysis, and archiving
- One validation saves time over validation of multiple instruments

Flexible Modular Design

BioProfile FLEX can be customized with 1–4 modules to consolidate up to 15 vital cell culture tests:

- Module 1: Gluc, Lac, Gln, Glu, NH_4^+ , pH, PCO_2 , PO_2 , Na^+ , K^+ , Ca^{++} , by Electrochemistry
- Module 2: Osmolality, by Freezing Point Depression
- Module 3: Cell Density/Cell Viability, by Digital Imaging
- Module 4: Optional Chemistries -- Phosphate by Photometry
Under Development

Gluc	Lac	Gln	Glu	NH_4^+	pH	PO_2	PCO_2	Na^+	K^+	Ca^{++}	CD	CV	Osm	PO_4
------	-----	-----	-----	-----------------	----	---------------	----------------	---------------	--------------	------------------	----	----	-----	---------------

nova
biomedical

IN THE U.S., call toll-free 800-458-5813 or 781-894-0800 • IN CANADA, call toll-free 800-263-5999
www.novabiomedical.com/gen0908.htm